


PRODUCT SPECIFICATION GUIDE

产品规格书

耳背式助听器

BTE Hearing aids

5124L

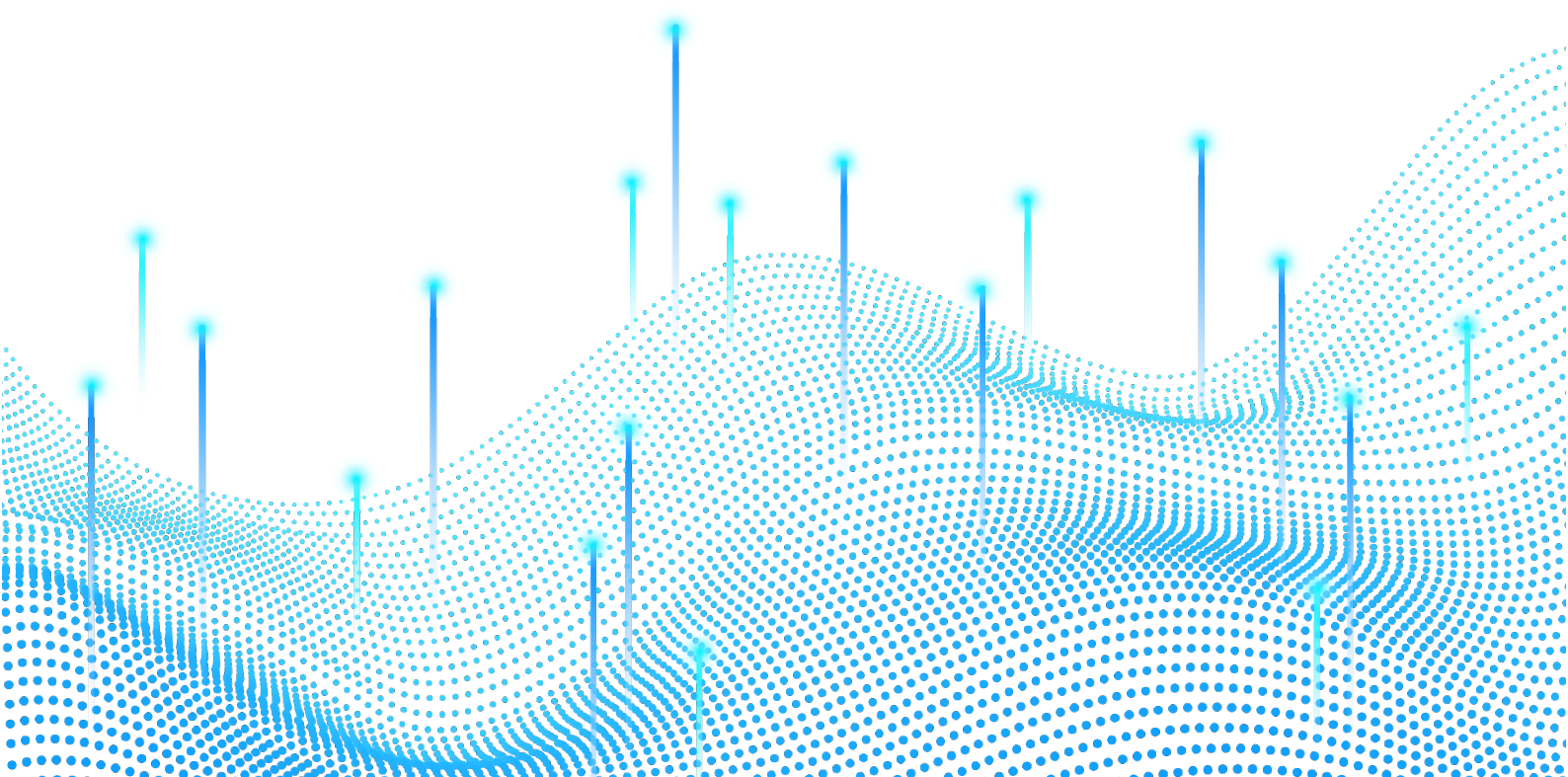

基本参数

BASIC PARAMETERS



双耳 Binaural

● 5124L

最大输出	Maximum Output	$\leq 114+3\text{dB SPL}$
平均输出	Average Output	$107\pm 4\text{dB SPL}$
峰值增益	Peak Gain	$40\pm 5\text{dB}$
平均增益	Average Gain	$35\pm 5\text{dB}$
频率范围	Frequency Range	$\leq 250\text{Hz} \sim \geq 5000\text{Hz}$
总谐波失真	Total Harmonic Distortion	$\leq 3\%$
等效输入噪声	Equivalent Input Noise	$\leq 20\text{dB}$
工作电流	Working Current	$\leq 3\text{mA}@3.7\text{V}$

WDRC通道数	Number of WDRC channels	16
电池类型	Battery Type	锂电池 Lithium Battery
麦克风类型	Microphone Type	硅麦 Silicon Microphone
喇叭类型	Speaker Type	动铁 Dynamic Iron Speaker
充电接口	Charging Port	TYPE- C

支持品牌LOGO 激光镭雕 Brand logo laser engraving is supported

注：默认黑白灰度镭雕，如需彩色丝印可咨询详谈

Note: Black-and-white grayscale laser engraving is provided by default. For color silk printing, please inquire for details.

- ▶ 提供效果图及打样定稿服务
Provide effect renderings and proofing finalization services.
- ▶ 根据LOGO样式，提供最优的OEM方案
Based on the logo design, we will offer the optimal OEM solution.



全新功能

NEW FEATURES

*自带消毒灭菌和干燥功能

Built-in sterilization, disinfection, and drying functions



紫灯亮起, 正在紫外线灭菌

The purple light is on, indicating ultraviolet sterilization in progress.




红灯亮起, 正在干燥储存助听器

The red light is on, indicating drying storage for the hearing aid.

增益调节


GAIN ADJUSTMENT

(1) 增益曲线 Gain Curve

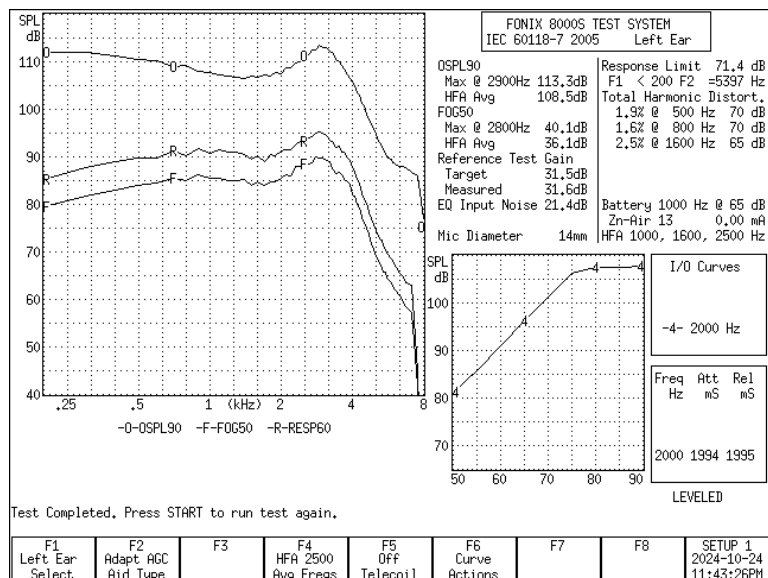

建议验配范围:

Recommended Fitting Range:

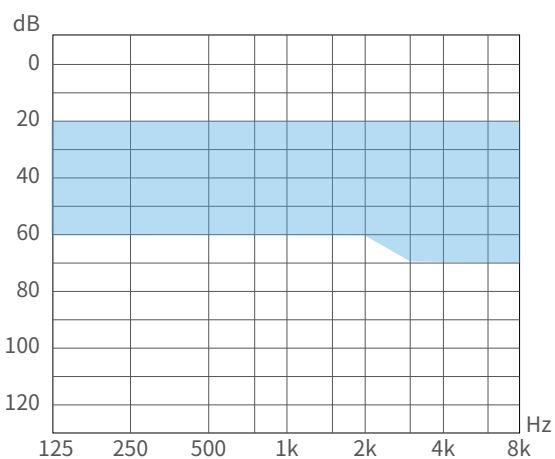

(2) 增益调节范围 Gain Adjustment Range

5124L: 10dB ~ 35dB

1档 Level 1	2档 Level 2	3档 Level 3	4档 Level 4	5档 Level 5	6档 Level 6	7档 Level 7	8档 Level 8	9档 Level 9
11dB	14dB	17dB	20dB	23dB	26dB	29dB	32dB	35dB

功能介绍

FUNCTION INTRODUCTION

(1) 蓝牙信息 Bluetooth Information

类型 Type: BLE 版本 Version: V5.3

(2) 操作说明 Instructions

	开机 Power on	关机 Power off	重启 Reboot manually
	充电盒中取出 Take out of the charging case	放回充电盒 Put the hearing aids back into the charging case	B键长按6秒 Press and hold the down button for 6 seconds
	9档音量 Level 9 Volume	切换蓝牙/助听 Switch hearing aid mode	切换助听模式 Switch hearing aid mode
	短按A/B Short press A to Volume up B to Volume down	按键B长按3秒 Press B and hold the button for 3 seconds	按键A长按3秒 Press A and hold the button for 3 seconds

(3) 续航 Battery life

 助听模式 Hearing mode 单次 Single: 15Hrs 综合 Comprehensive: 77Hrs	
 蓝牙模式 Bluetooth mode 单次 Single: 4Hrs 综合 Comprehensive: 20Hrs	
双耳 Binaural 充电盒电池容量 Charging case battery capacity: 400mAh	
助听器电池容量 Hearing aid battery capacity: 30mAh	
耳机充电时间 Charging time (Earphones): 2Hrs 充电仓充电时间 Charging time (Case): 2Hrs	充电仓给助听器可供充电次数: 4Times Number of charging cycles the case supports

(4) APP验配 APP Fitting

- FIG6 验配公式 FIG6 Fitting Formula
- 4个模式操控调节 (通用、降噪、耳鸣掩蔽、蓝牙) 4 mode control adjustment (General, Noisy, Tinnitus Masking, Bluetooth)
- 音量调节: 左耳调节、右耳调节、APP内左右耳同步调节 Volume adjustment: left ear adjustment, right ear adjustment, In-app left-right ear synchronization adjustment.
- 支持电量百分比查看 Support battery percentage viewing (ForSound)
- 用户可以调整大声、小声、低音和高音的音调 (32 频段) The user can adjust the tone of louder, softer, bass, and treble. (32 band)

FIG 6 国际验配公式

采用门店验配师大量使用的FIG 6验配公式进行目标曲线的验配,达到线下专业级要求

Download on the App Store

GET IT ON Google Play


FIG6 International Fitting Formula:
Achieves professional-grade offline requirements by calibrating target curves with the FIG6 formula widely adopted by in-store audiologists.

(5) 芯片介绍 Chip Introduction

Onsemi 进口专业级助听器芯片 Imported Professional-grade Hearing Aid Chips

• 自适应宽动态范围压缩算法 Adaptive Wide Dynamic Range Compression	
• 噪声抑制算法 Environmental Noise Reduction	• 声反馈抑制算法 Digital Feedback Canceller
• 瞬态强声阻断 Instantaneous Strong Sound Interruption	• 风噪声抑制 Wind Noise Reduction
• 自适应方向性麦克风 Adaptive Directional Microphone	• 音量学习 Volume Learning
• 耳鸣掩蔽 Tinnitus Masking	• 数据记录 Data Logging
• 自适应级别管理 Adaptive Level Management	• 实时分析演示 Dynamic Display

(6) 支持PC验配 PC Fitting

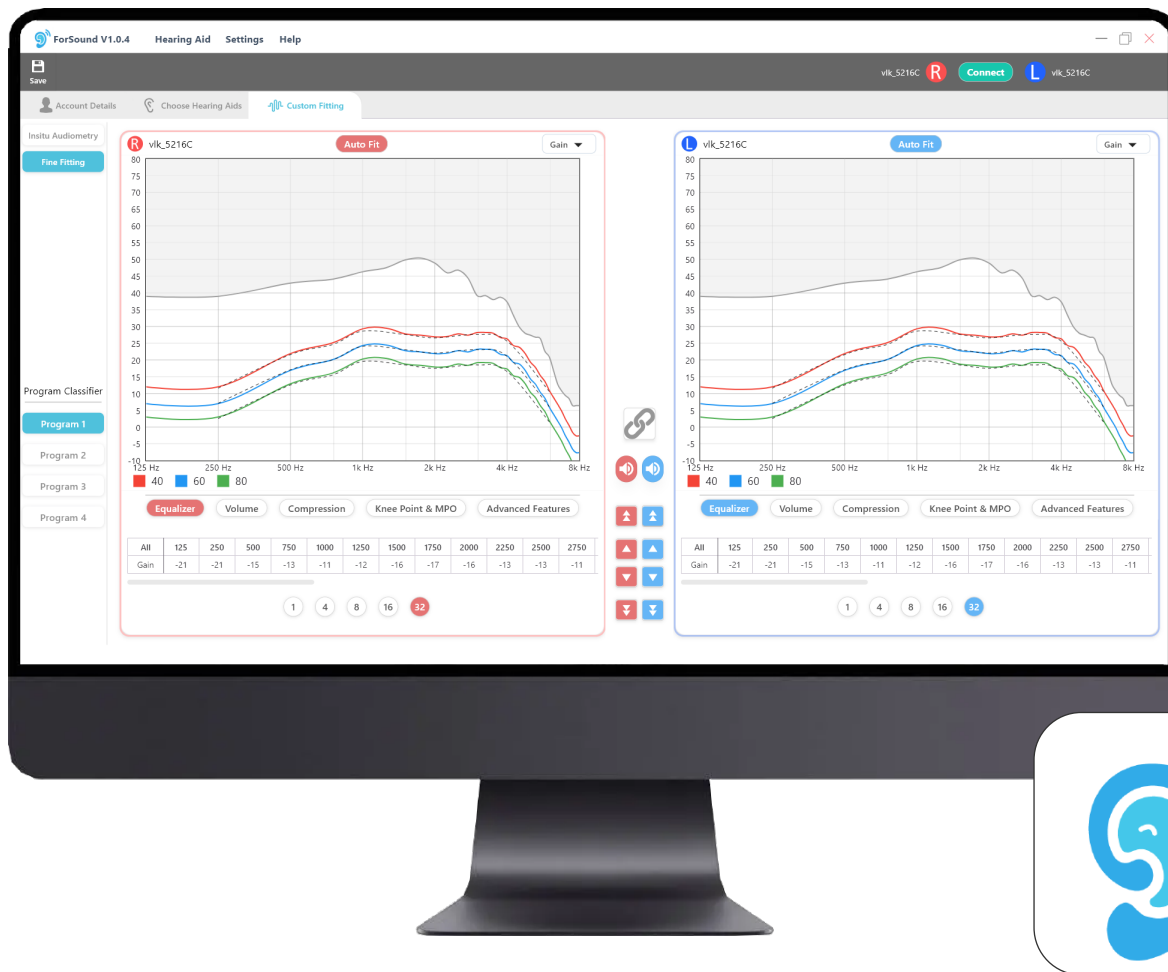

- 支持HiPro、Communication Accelerator Adaptor、Builted Ble、Ble Dongle 连接
Supports HiPro、Communication Accelerator Adaptor、Builted Ble、Ble Dongle connection
- FIG6 验配公式 FIG6 Fitting Formula
- 根据用户需求自定义验配 Custom fitting based on user needs
- 音量、压缩比、拐点和 MPO均可调整方式同均衡器调整
Volume, compression ratio, knee, and MPO are adjustable in the same manner as equalizer settings
- 可保存用户信息和调整数据以便随访和后续调整
Stores user information and adjustment data for follow-up and subsequent modifications

产品外观

PRODUCT APPEARANCE



● 珠光黑+枪灰

Pearl Black + Gunmetal Gray



● 浅金+枪灰

Soft Gold + Gunmetal Gray

配件清单

ACCESSORIES LIST



助听器耳机
Hearing aids




助听器充电仓
Charging Case



L

S

双层开放式耳帽
Open Double-Layer Dome

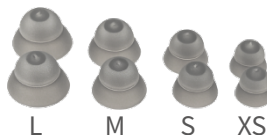


L

M

S

XS

双层耳帽
Double-Layer Dome

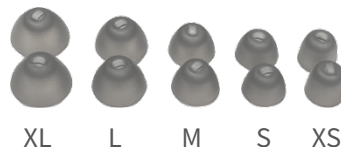

XL

L

M

S

XS

单层耳帽
Single-Layer Dome



USB-TypeC 数据线 (1m)
USB-C Cable (1m)




耵聍帽*8
Wax Guard*8



清洁刷
Cleaning Brush



说明书
User Manual

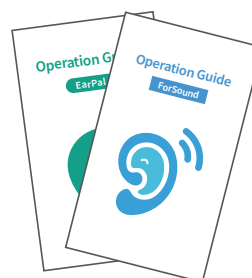

APP操作指南
App Operation Guide




提示卡
Tips card

尺寸与重量

SIZE AND WEIGHT

双耳 Binaural

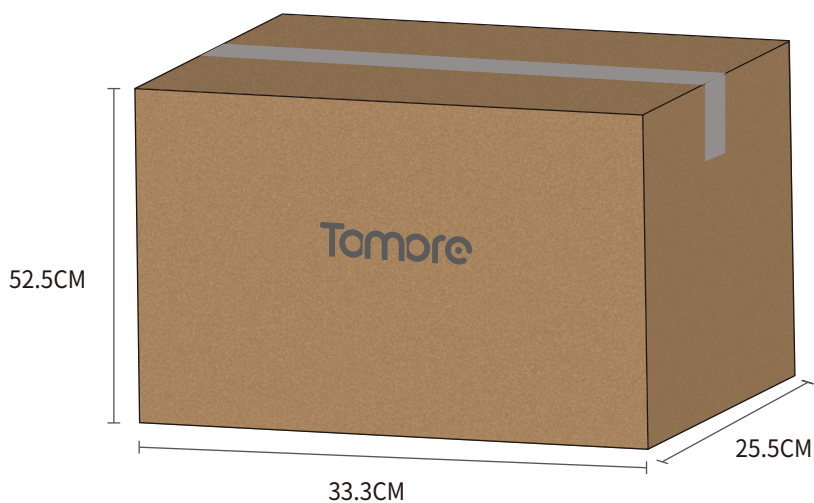

	长 Length	宽 Width	高 Height	重量 Weight	50台/箱 Pcs/ctn
助听器 Hearing aid	12.3mm	7.9mm	28.1mm	3.2g/Single	
充电盒 Chging case	84mm	31mm	58.8mm	69.5g	
包装盒 Packaging box	160mm	100mm	46mm	290g	
外箱 Outer Box	52.5cm	33.3cm	25.5cm	15.25kg	